


# Emergency Caller Products

GETS  
HELP  
EAST


(T) 800-227-2474 (F) 800-762-7613

## CC-900 Wireless Enunciator with Pendant QUICK START GUIDE


Note\* For use with  
Emergency Caller Products Wireless  
900MHz Supervised Transmitters Only!

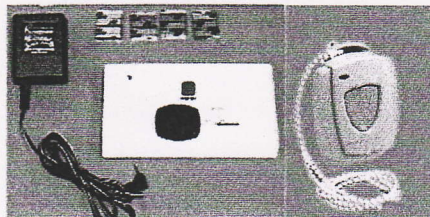


- 1 Components:
  1. Enunciator (Receiver Horn-Strobe)
  2. Power Adapter
  3. Four (4) Velcro tabs
  4. Pendant with necklace

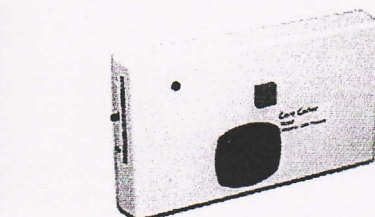

- 2 The Care Caller system is designed to assist the caregiver in monitoring patient activity. It is selectable for LOW or HIGH horn sound and comes with a Flashing Strobe. Selectable Off position when not in use. Unit has Aux. out (2.5mm jack, normally open) for Automatic Voice Dialer other connections.

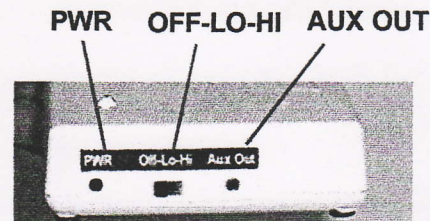


- 3 Keep in mind the Enunciator has two sound levels- low and high. Also remember the transmitter range will be dependent to some degree on the structure of the residence. Also, you will need to place the unit within reach of an electrical outlet in order to plug in the adapter.

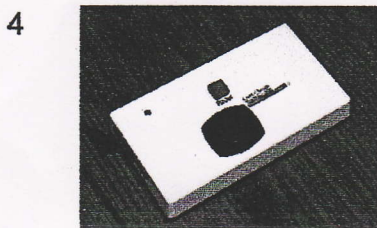

The Enunciator may be placed at Nurse Station, desk, on counter, nightstand, wall mounted or outside (covered). The caregiver should also consider where he or she might frequently be to be within sound of the enunciator when it alarms.

**Important!** Be sure to test your system often, with transmitter at locations where patient may be.

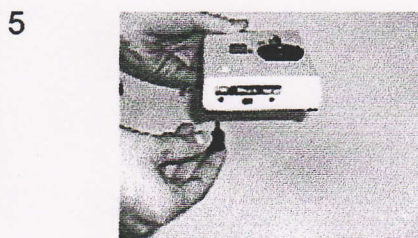

Now that we have decided where to place the enunciator let's go on to the next step. Begin by plugging the adapter into the top of the Enunciator by PWR as shown.


Next, plug the power adapter into the wall socket.

The pendant has been pre-programmed to work with the enunciator, so no action is needed for it at this point.



Now you are ready to test the Enunciator. Move the switch on the side of the unit to either low or high, for the sound output.

**Test the system.** Activate the transmitter and listen for the sound. If you do not hear the Enunciator sound, check to be sure the switch is in Lo or Hi position and that the adapter is plugged in the wall. When your Enunciator sounds you have successfully setup and tested your system.

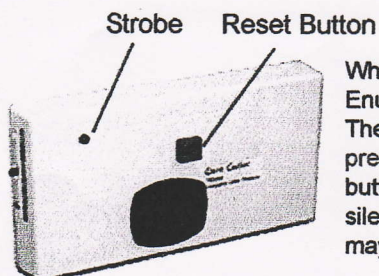

When activated the transmitter will send a signal to the Enunciator and the alarm sounds and the strobe flashes. The Enunciator continues to sound until reset button is pressed by Caretaker, assuring Help Has Arrived! Hold the button down for at least ten (10) seconds. The sound is then silenced, until transmitter is activated again. The Enunciator may be turned off when operation not desired.

Congratulations! Your system is all set up. Now you can feel Safe, Secure and at Ease whenever patients, family members, or friends are confined to bed, home or facility and need fast assistance. Care Caller sends the signal to call for help up to 1000 feet indoors (depending on structure, metal obstacles etc.). Be sure to test your system often as range may be effected by conditions.

**Note: The enunciator emits beeping sound whenever Emergency Caller Transmitter battery is getting low.**

See Reverse side for a list of other ECP Products.